

10x Passivation solution

This solution is intended for use as a passivation solution on a **switch**SENSE® chip. It contains a thiol-reactive compound for the formation of a self-assembled-monolayer.

Product Description

Order Number **SOL-PAS-1-5**

Material	Cap	Amount	Storage	Comment
10x Passivation solution; 0.2 µm sterile filtered	yellow	5 x 1 mL	-20°C	

For *in vitro* use only.

This product has a limited stability, please see expiry date on label.

Preparation

Dilute 10x Passivation solution with deionized 0.2 µm filtered water by mixing the complete solution (1 mL) with 9 mL ultrapure water.

The diluted solution should be stored at 2-8°C and can be used for 2 weeks.

Contact

Dynamic Biosensors GmbH

Perchtinger Str. 8/10
81379 Munich
Germany

Phone: +49 89 89 74 544 0

Dynamic Biosensors Inc.

300 Trade Center, Suite 1400
Woburn, MA 01801
USA

Phone: +1 781 404 6126

Order Information


Phone: +49 89 89 74 544 0

Email: order@dynamic-biosensors.com

Technical Support

Phone: +49 89 89 74 544 66

Email: support@dynamic-biosensors.com

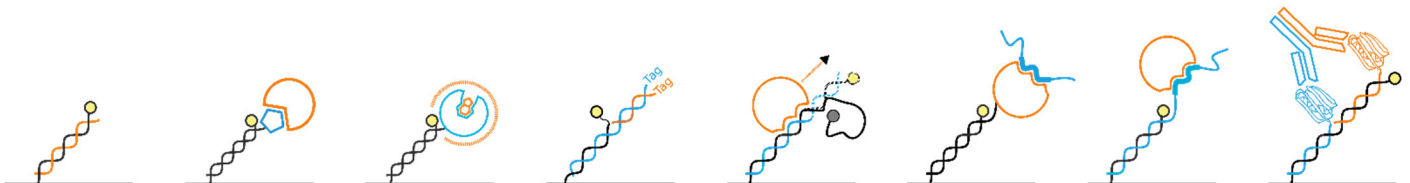

switchSENSE® is a proprietary measurement technology by Dynamic Biosensors GmbH. Instruments and biochips are engineered and manufactured in Germany.

©2023 Dynamic Biosensors GmbH | Dynamic Biosensors Inc. All rights reserved.